


Ghost Zero Column Application : Clopidogrel bisulfate [USP][유연물질]


System Connection

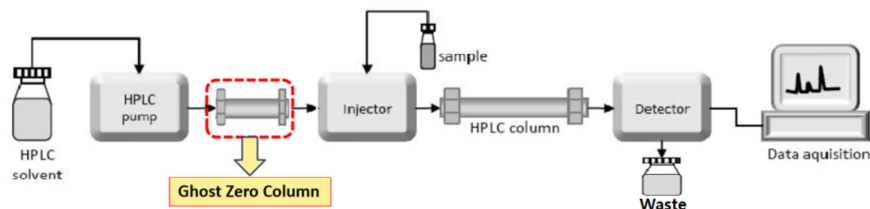

Application and Results

Mobile Phase A	0.96 g/L sodium 1-pentanesulfonate. Adjust with phosphoric acid to a pH of 2.5
Mobile Phase B	Acetonitrile
Mobile Phase C	Methanol
Flow Rate	1 mL/min
Temperature	30 °C
Detector	UV 220 nm
Injection Volume	10 µL
Sample	Acetonitrile : Buffer = 60 : 40 Solution
Column	INNO Column C18 5µm 4.6 * 150mm

Time (min)	Mobile Phase A (%)	Mobile Phase B (%)	Mobile Phase C (%)
0	85	10	5
3	85	10	5
48	30	65	5
68	30	65	5
78	85	10	5

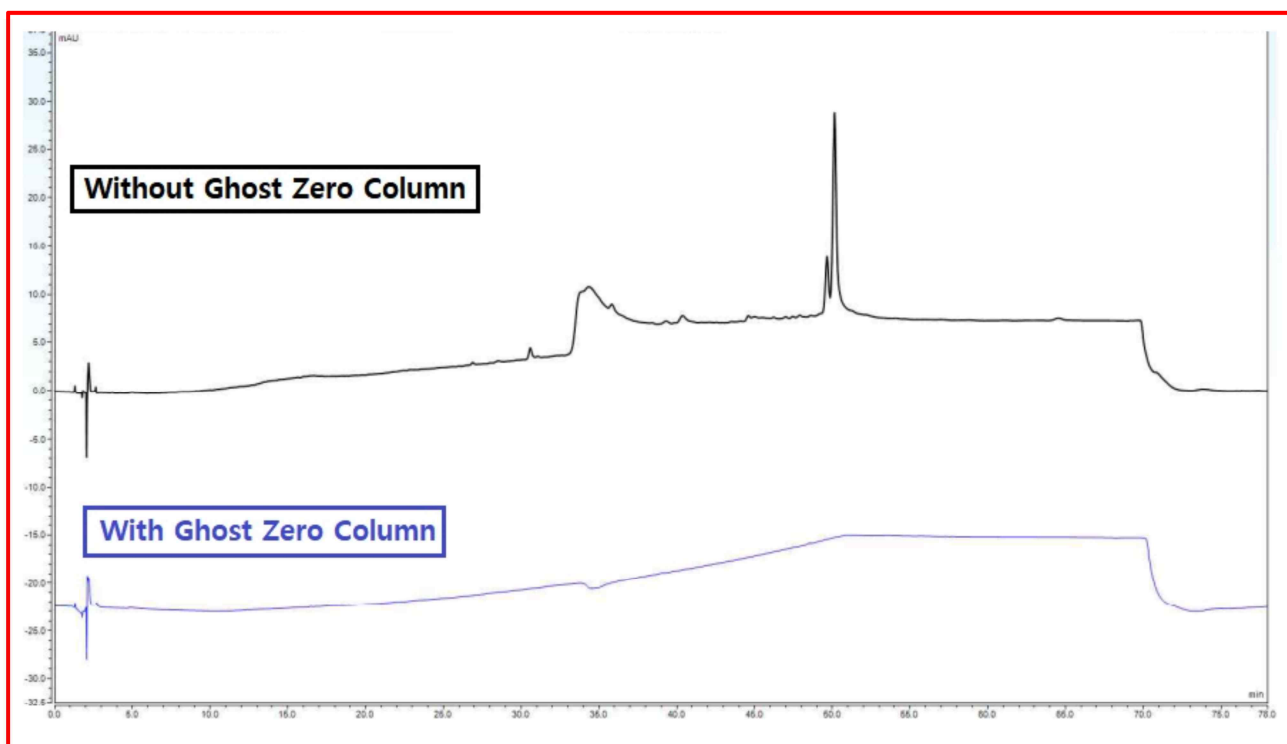

제품 정보

제품명	Ghost Zero Column	규격	4.6 * 50mm	제품 번호	GZ04605
-----	-------------------	----	------------	-------	---------