


HPLC DATA SHEET

Amino acids (15 kinds)

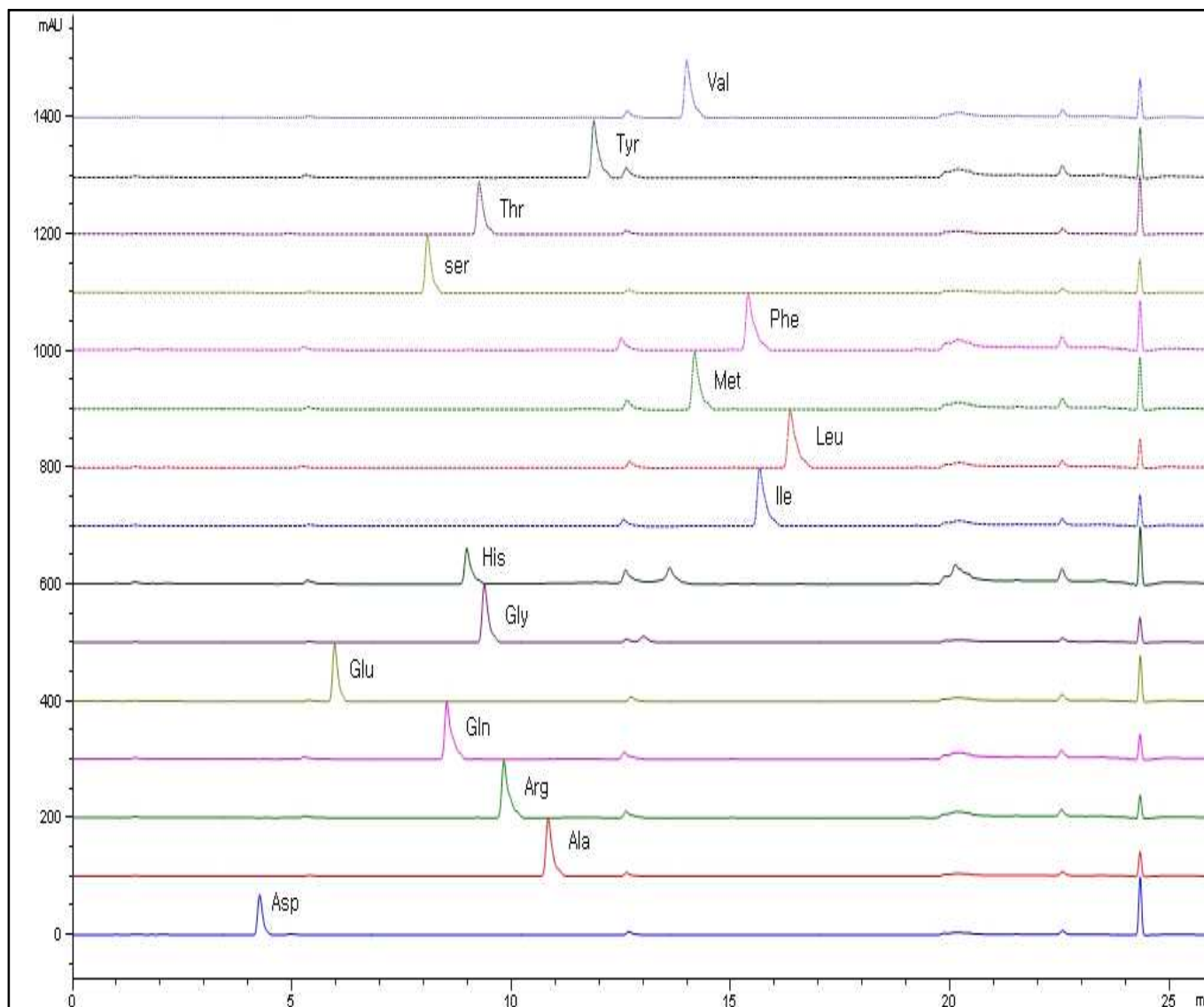

HPLC condition

Column INNO Column C18 (5 μ m, 4.6mm I.d. \times 150mm)
Flow rate 1.5 mL/min
Temperature 40 $^{\circ}$ C
Detector UV (338nm)
Mobile Phase (A) : 40mM Sodium phosphate(di-basic) 0.1% phosphoric acid
Mobile Phase (B) : ACN / Methanol / DW = 45 / 45 / 10
Mobile Phase (B) 0%(0min)-57%(18.1min)-100%(18.6min)-100%(22.3min)-0%(23.2min)-0%(26min)
Gradient