


HPLC DATA SHEET

Vitamins – Water soluble

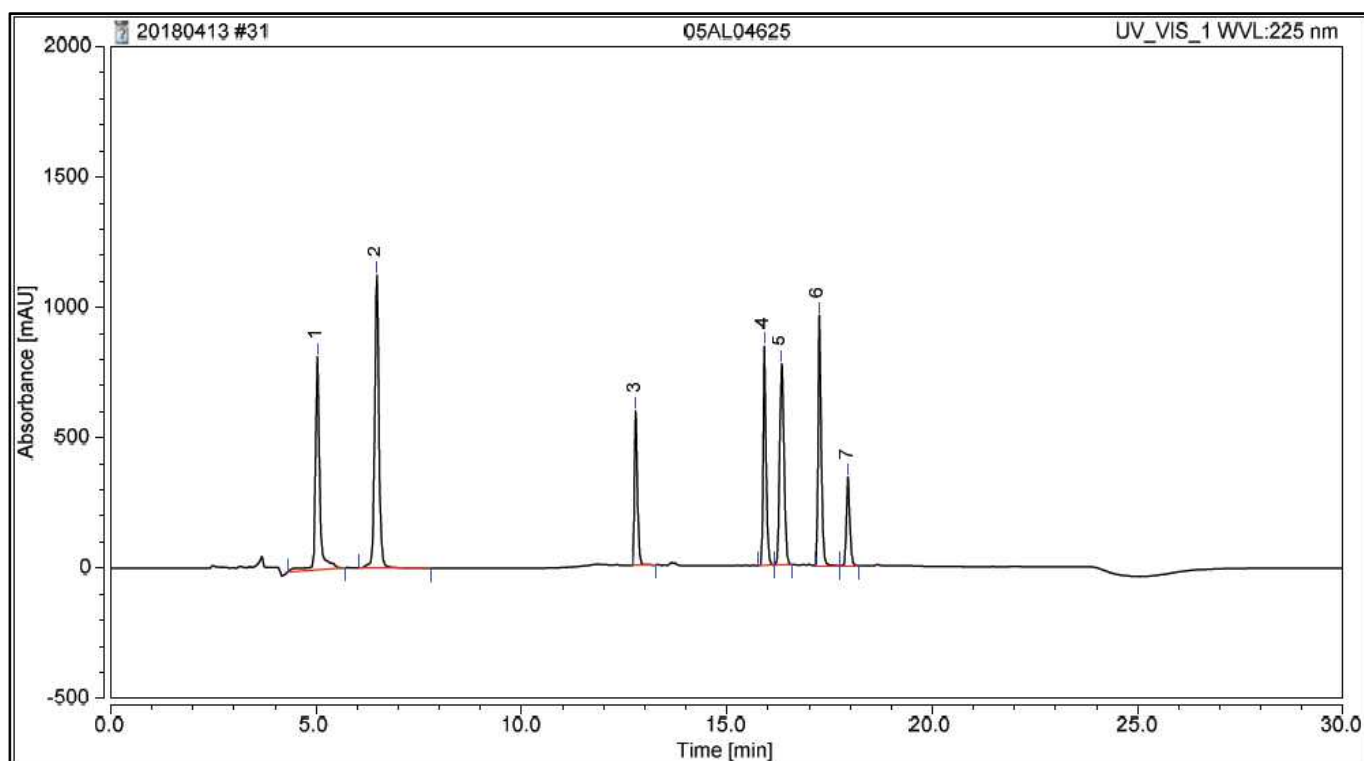

1. Ascorbic acid
2. Nicotinamide
3. Pyridoxine HCl
4. Folic acid
5. Cyanocobalamin(B12)
6. Riboflavin
7. Biotin

HPLC condition

Column Aegispak C18-L (5um, 4.6mm I.d. × 250mm)

Mobile Phase (A) : 0.05% TFA (Water)

Mobile Phase (B) : 0.05% TFA (ACN)

(B) 1%(0min)-1%(5min)-30%(15min)-30%(20min)-1%(20.1min)-1%(25min) Gradient

Flow rate 0.8mL/min

Temperature 35°C

Detector UV (225nm)