


HPLC DATA SHEET

Vitamins – Water soluble

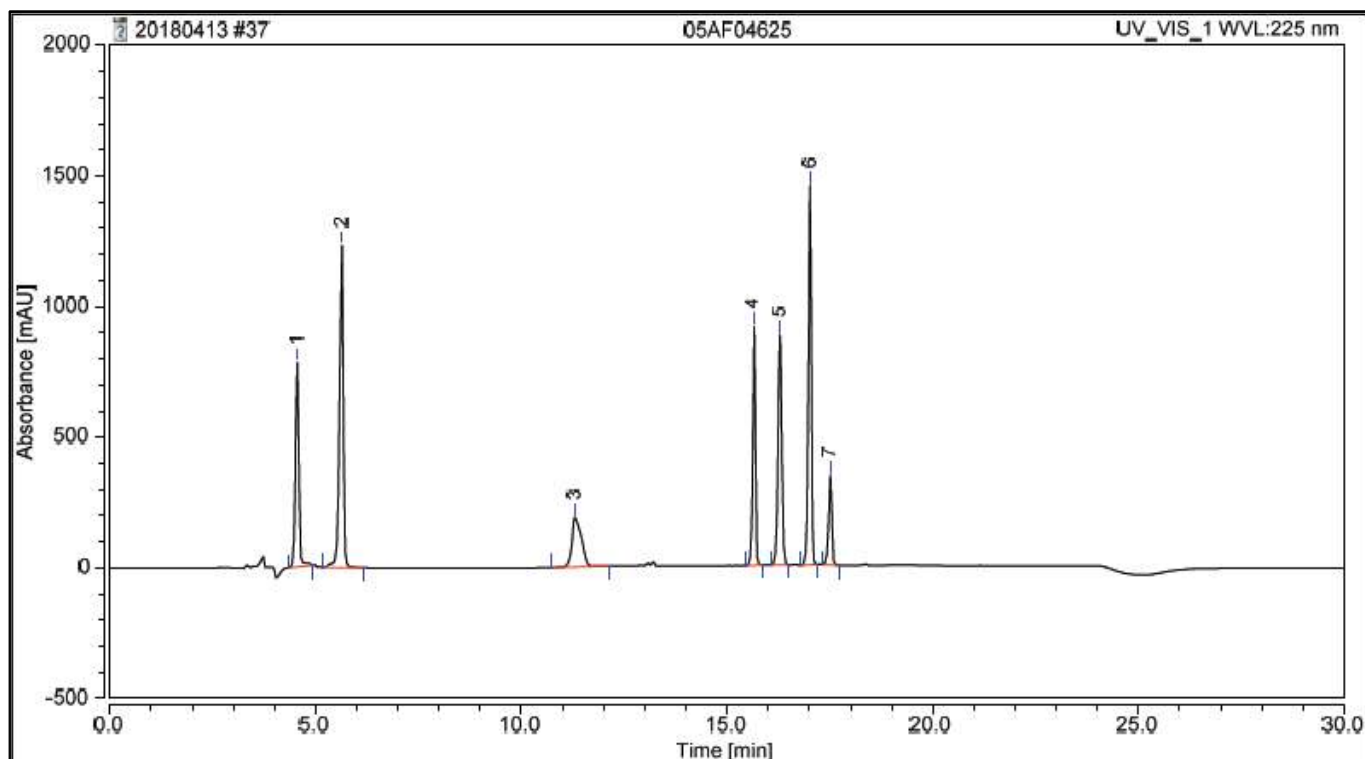

1. Ascorbic acid 2. Nicotinamide 3. Pyridoxine HCl 4. Folic acid
5. Cyanocobalamin(B12) 6. Riboflavin 7. Biotin

HPLC condition

Column	Aegispak C18-F (5um, 4.6mm I.d. × 250mm)
Mobile Phase	Mobile Phase (A) : 0.05% TFA (Water)
	Mobile Phase (B) : 0.05% TFA (ACN)
	(B) 1%(0min)–1%(5min)–30%(15min)–30%(20min)–1%(20.1min)–1%(25min) Gradient
Flow rate	0.8mL/min
Temperature	35℃
Detector	UV (225nm)