


HPLC DATA SHEET

Cyanocobalamine (Multivitamin Tablets)

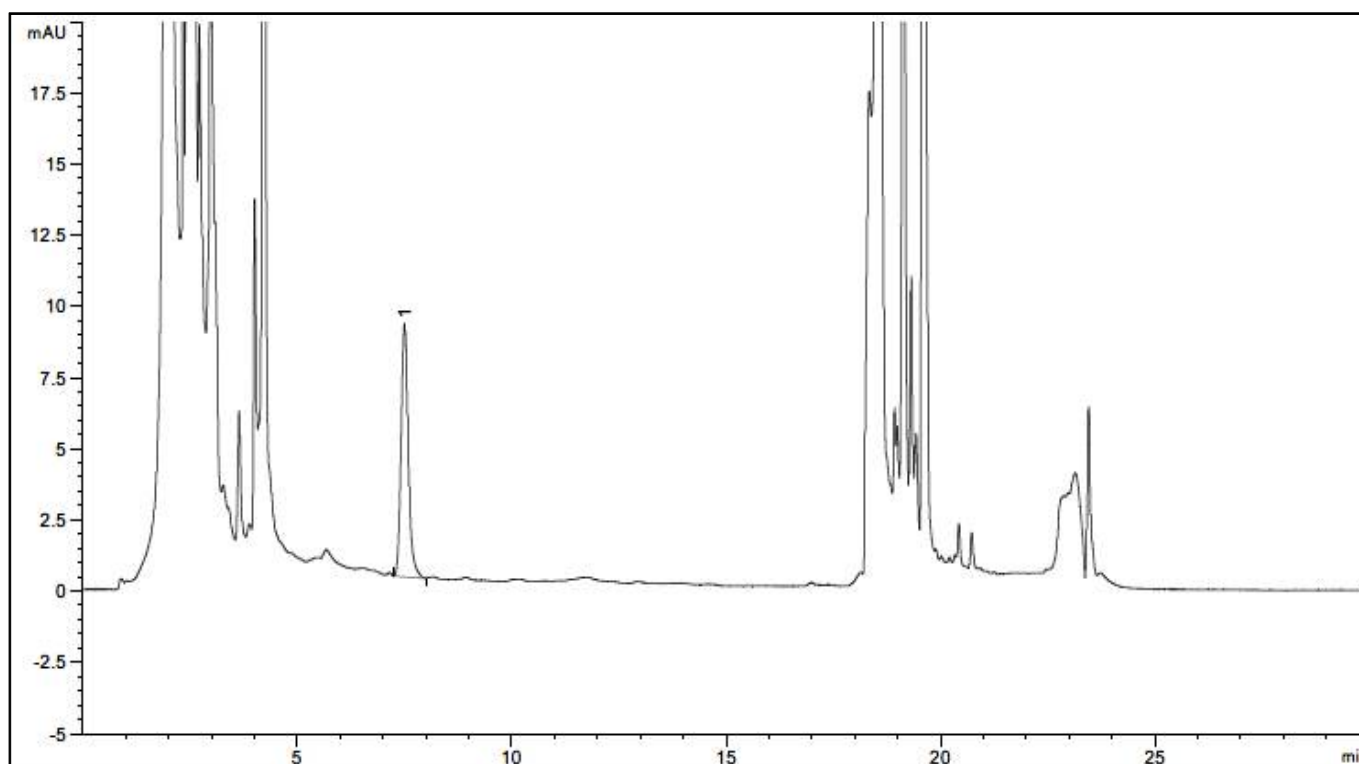

1. Vitamin B12 (Cyanocobalamine) : Multivitamin tablet

HPLC condition

Column INNO Column C18 (3 μ m, 4.6mm I.d. \times 150mm)

(A) Water

(B) Methanol

Mobile Phase

Time(min)	(A, %)	(B, %)
0.01	75	25
15.0	75	25
15.1	0	100
20.0	0	100
20.1	75	25
30.0	75	25

Flow rate 0.8 mL/min

Temperature 35 $^{\circ}$ C

Detector UV (360nm)