


HPLC DATA SHEET

Tocopherol Acetate (KP, Assay)

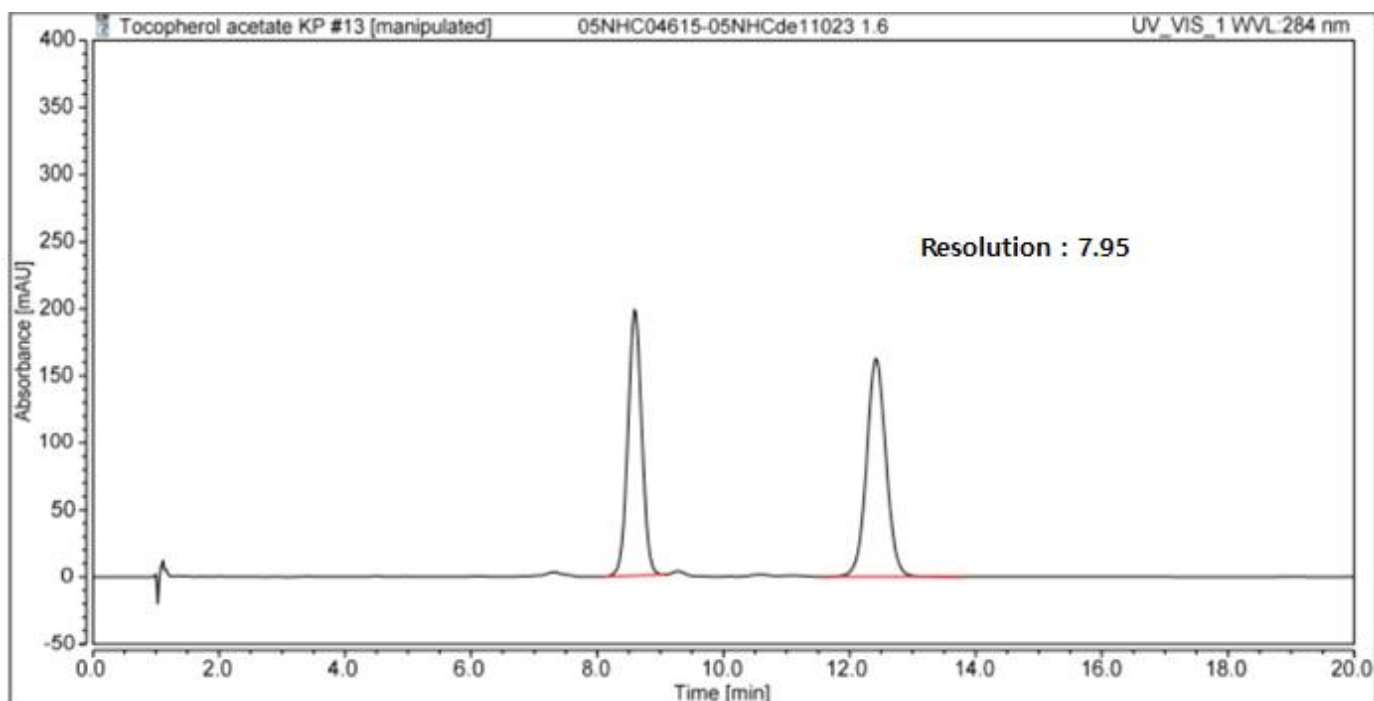

1. Tocopherol
2. Tocopherol Acetate

HPLC condition

Column INNO Column C18-HC 5um, 4.6mm I.d. × 150mm (P/N = 05NHC04615)

Mobile Phase Methanol : Water (49 : 1)

Flow rate Adjust so that the retention time of Tocopherol Acetate is about 12 minutes
(About 1.6 mL/min)

Temperature 35 °C

Detector UV (284 nm)

Injection Volume 20 uL