


HPLC DATA SHEET

Cilostazol (KP, JP, Related Substances)

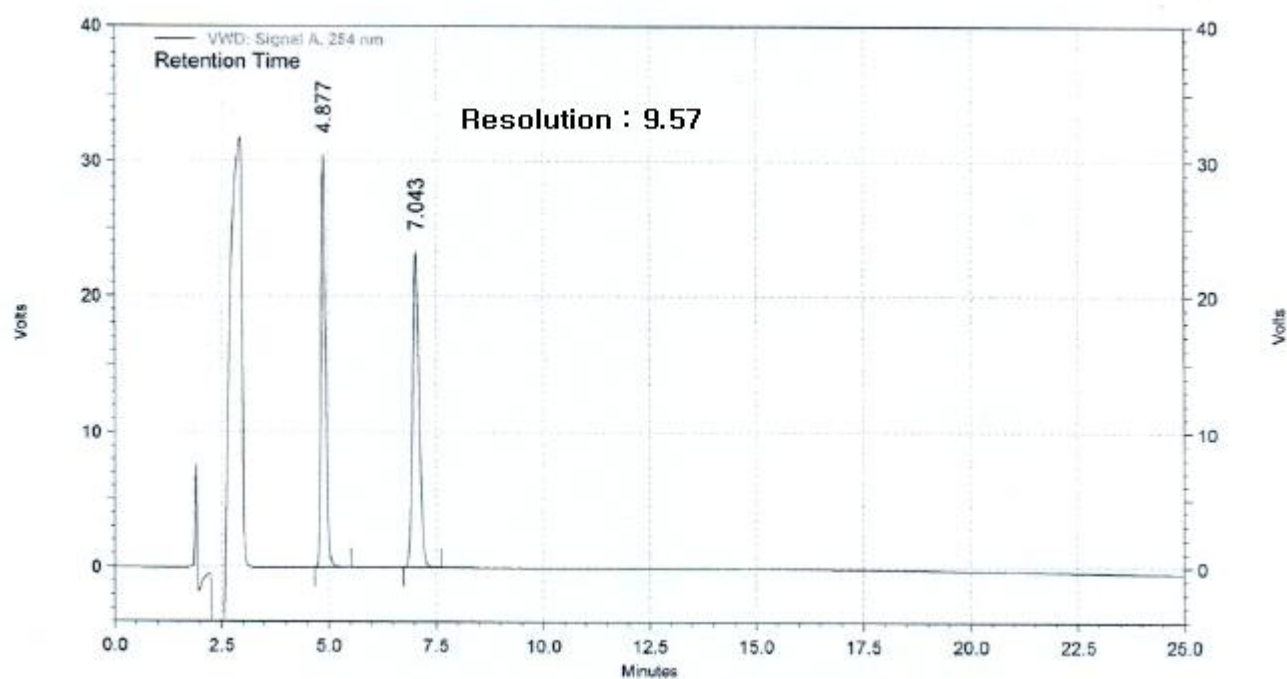

1. 3,4-dihydro-6-hydroxy-2(1H)-quinolinone 2. Cilostazol

HPLC condition

Column INNO Column Silica (5 μ m, 4.6mm I.d. $\times$ 150mm)

Mobile Phase Hexane : Ethyl acetate : Methanol = 10 : 9 : 1

Flow rate 1.0~1.1mL/min(Adjust so that the retention time of cilostazol is about 7 minutes)

Temperature 25 $^{\circ}$ C

Detector UV (254nm)