


HPLC DATA SHEET

Flecainide Acetate Tablets (USP)

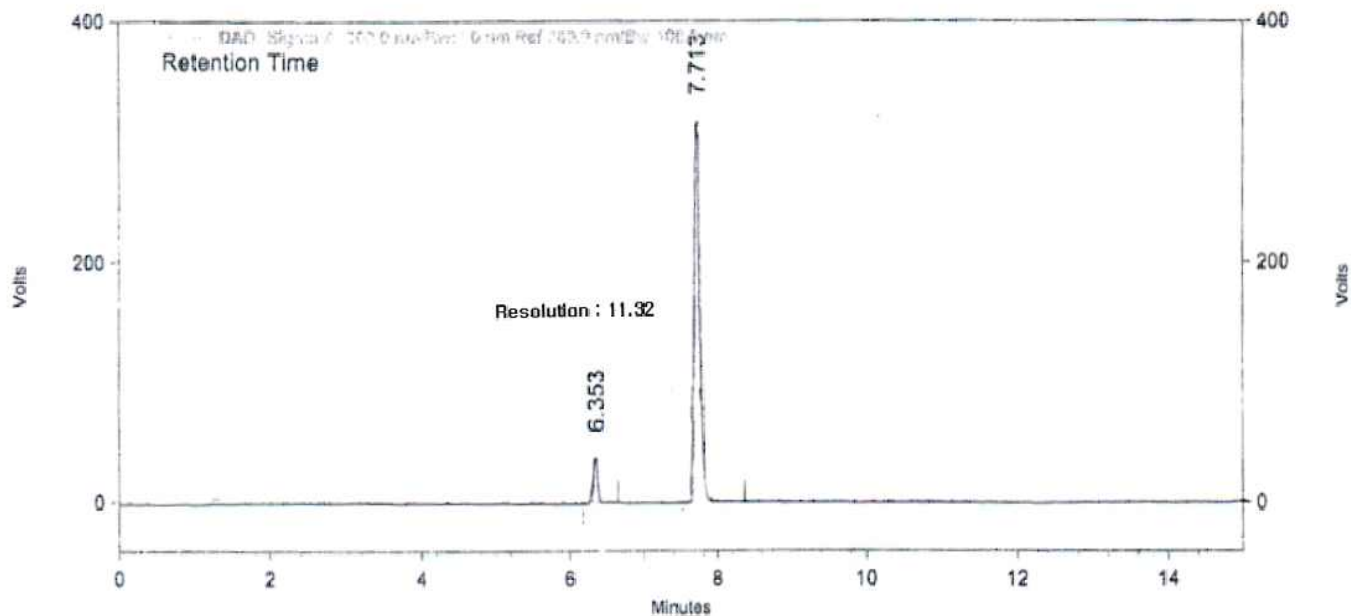

1. Flecainide Acid 2. Flecainide Acetate

HPLC condition			
Column	INNO Column C18 (3um, 3.0mm I.d. × 150mm)		
Mobile Phase	Mobile Phase A : 0.1M ammonium acetate (Adjust with acetic acid to a pH of 5.4)		
	Mobile Phase B : Acetonitrile		
	Time	Mobile Phase A	Mobile Phase B
	0	75	25
	7.0	45	55
	12.0	45	55
12.1	75	25	
15.0	75	25	
Flow rate	0.5mL/min		
Temperature	25°C		
Detector	UV (300nm)		