


HPLC DATA SHEET

Acyclovir Injection

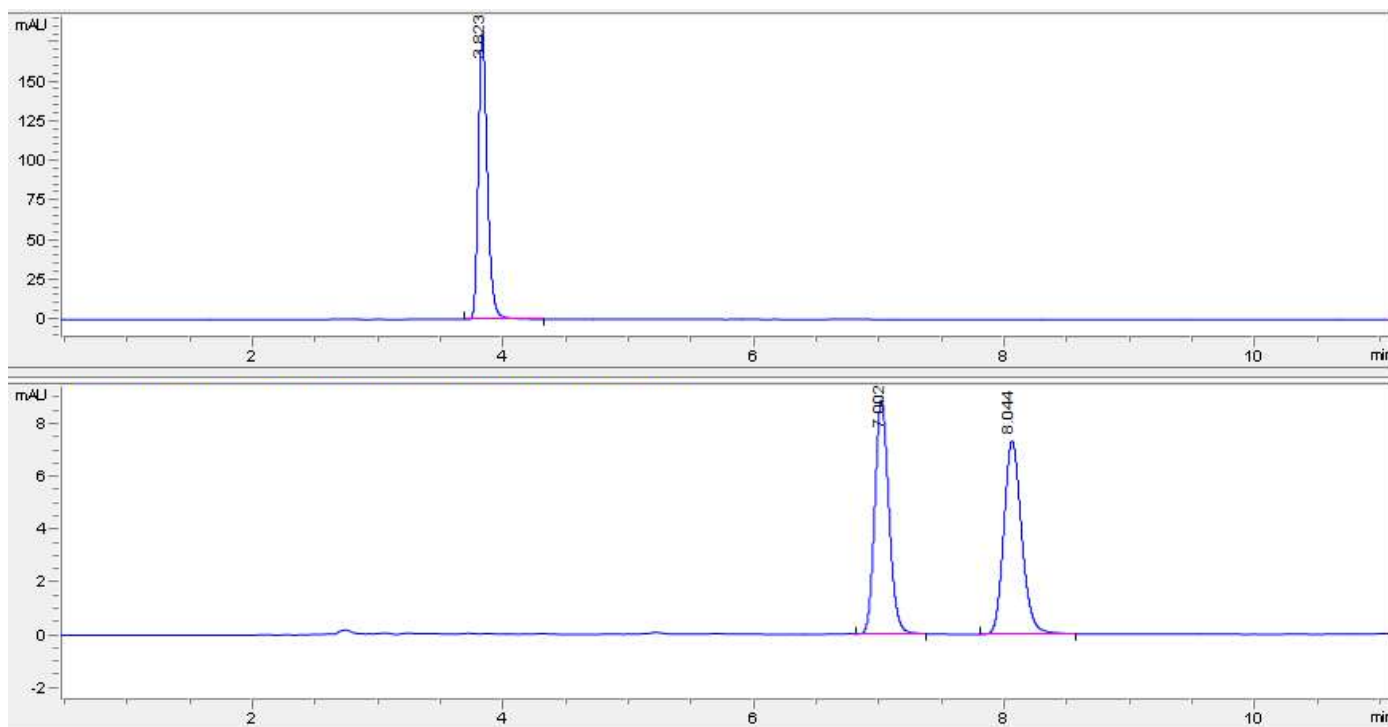

1. Guanine 2. Purine 3. Acyclovir

HPLC condition

Column INNO Column C18-SB (5um, 4.6mm I.d. × 250mm)

Mobile Phase A : 0.17M Acetic acid : Methanol = (125 : 8)

Mobile Phase B : Methanol

Mobile Phase	Time	Mobile Phase A	Mobile Phase B
	0	100	0
	15	100	0
	45	65	35
	46	100	0
	56	100	0

Flow rate 1.0mL/min

Temperature 25°C

Detector UV (254nm)