


HPLC DATA SHEET

Mupirocin (Assay)

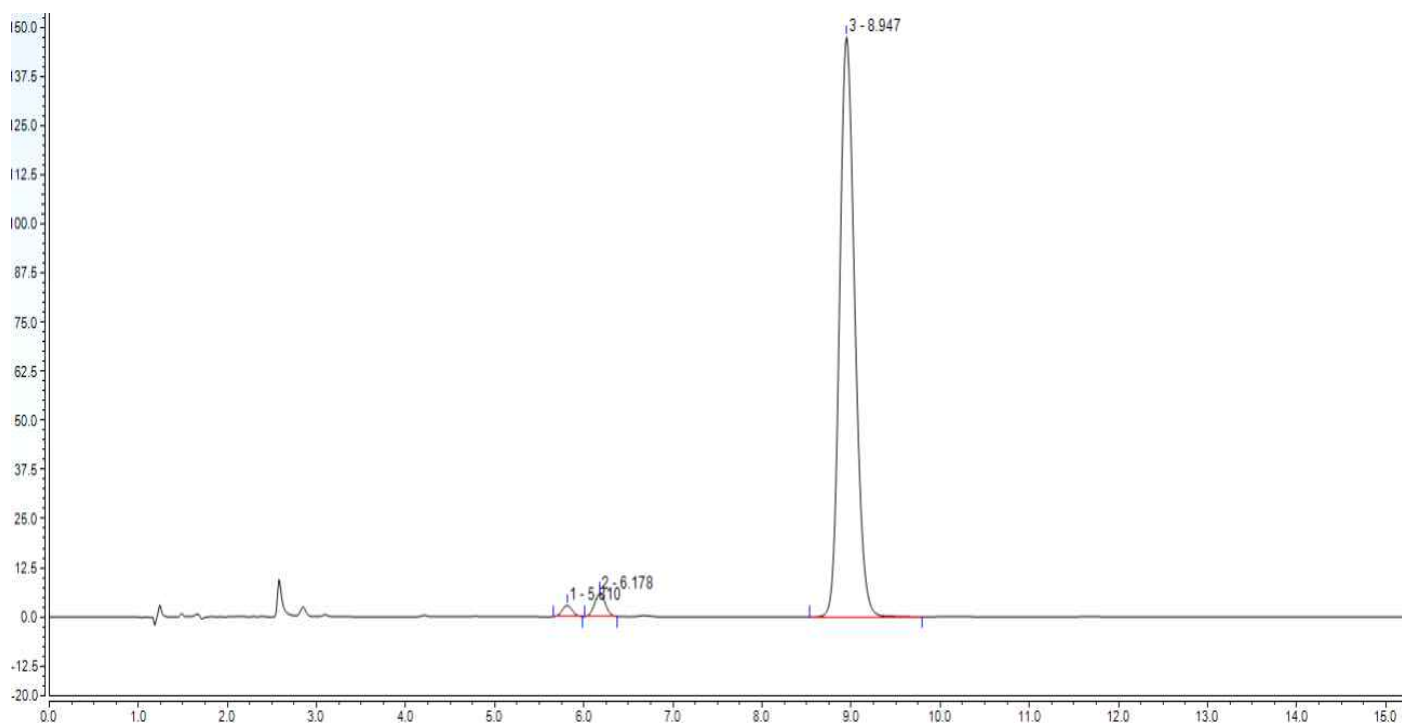

1. Mupirocin hydrolysate 2. Mupirocin hydrolysate 3. Mupirocin

HPLC condition

Column INNO Column C18-SB (5um, 4.6mm I.d. × 250mm)

Mobile phase : Solution A(pH6.3) : Acetonitrile (75 : 25)

Mobile Phase *Solution A(pH6.3) : 0.05mol/L Sodium dihydrogen phosphate monohydrate
(adjust the pH to 6.3 with 10 mol/L Sodium hydroxide)

Flow rate 2.0mL/min

Temperature 25°C

Injection Vol. 20uL

Detector UV (229nm)