


HPLC DATA SHEET

Cyclobenzaprine Hydrochloride Tablets (USP, Assay)

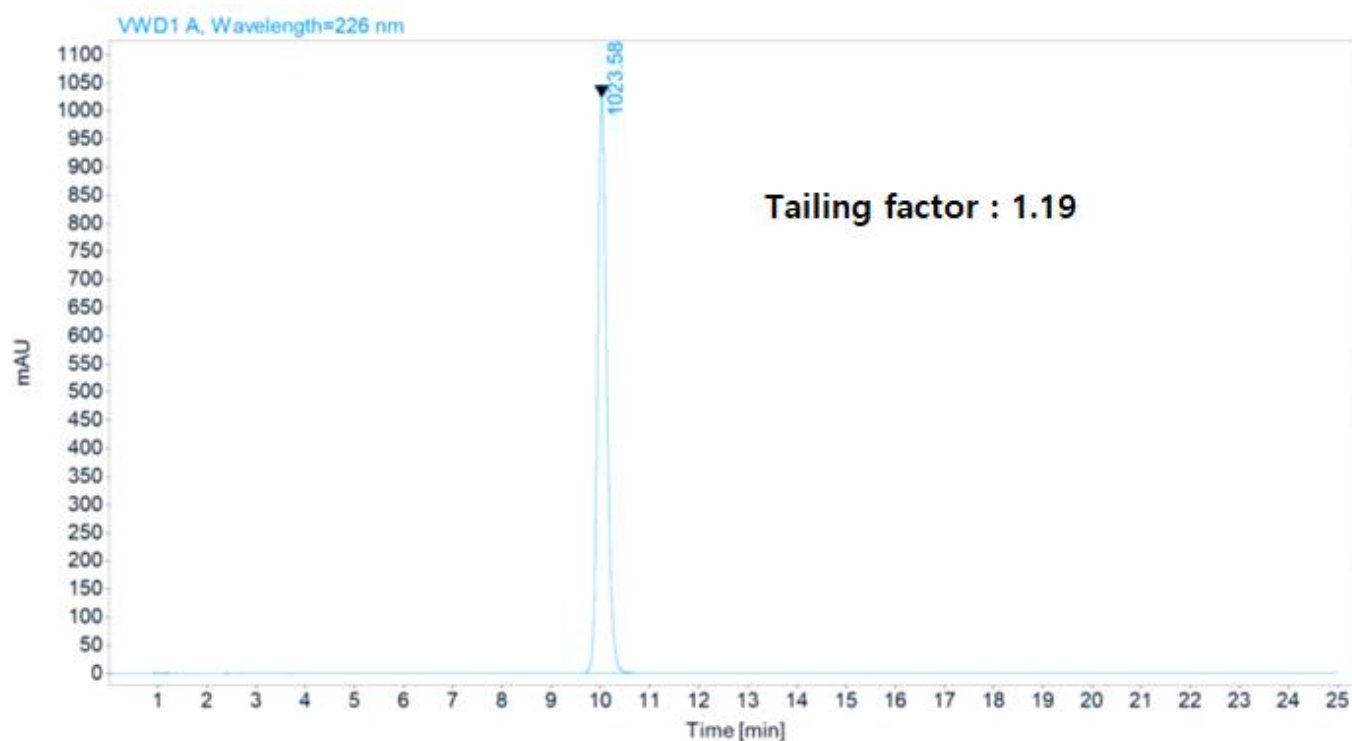

1. Cyclobenzaprine Hydrochloride

HPLC condition

Column	INNO Column C8-SB (5 μ m, 4.6mm I.d. \times 250mm)
Mobile Phase	Methanol : Buffer = (65 : 35) *Buffer : 11.4g/L of ammonium acetate in water. Adjust with ammonium hydroxide to a pH of 7.2.
Flow rate	1 mL/min
Temperature	30 $^{\circ}$ C
Injection Vol.	10 μ L
Detector	UV (226 nm)