


HPLC DATA SHEET

Tyrothricin (Related substances)

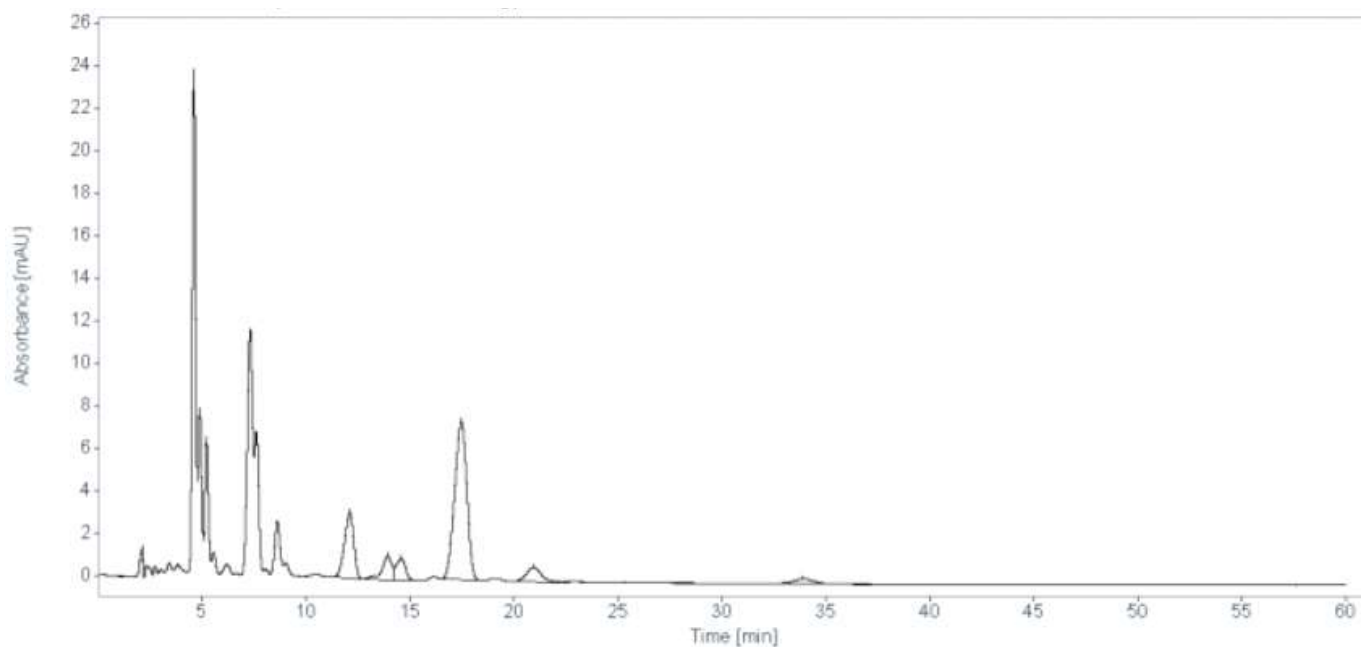

Sample = Tyrothricin

HPLC condition

Column INNO Column C18-HC 5 μ m, 4.6mm I.d. \times 250mm (P/N 05NHC04625)

Mobile Phase Buffer : Methanol = 75 : 25
* Buffer : Dissolve 0.79g of ammonium sulfate in 1L of water

Flow rate 1.2 mL/min

Temperature 60 $^{\circ}$ C

Detector UV (280 nm)

Injection volume 25 μ L