


# HPLC DATA SHEET

Copovidone (USP, 1-vinyl-2-pyrrolidone)

## Copovidone 1-Vinyl-2-Pyrrolidone USP

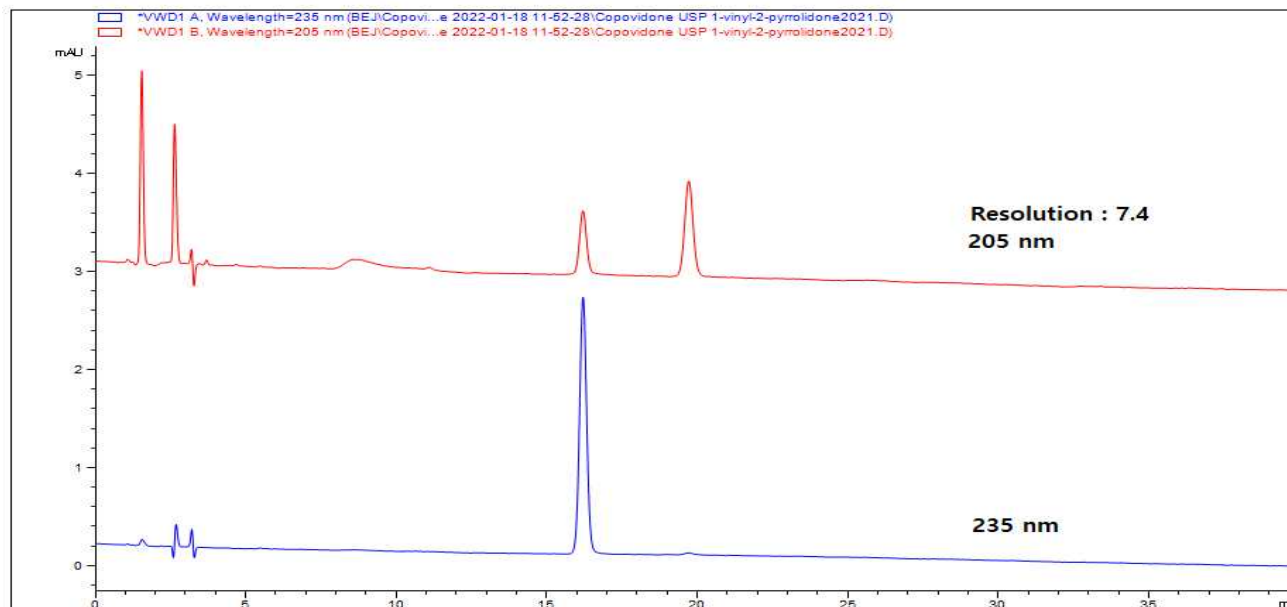

1. 1-vinyl-2-pyrrolidone      2. vinyl acetate

### HPLC condition

Column                    Aegispak C18-L 5 $\mu$ m, 4.0mm I.d.  $\times$  33mm (P/N 05AL040033)  
                                 Aegispak C18-L 5 $\mu$ m, 4.0mm I.d.  $\times$  250mm (P/N 05AL04025)

Mobile Phase            DW : Acetonitrile = 23 : 2

Flow rate                1.0 mL/min

Temperature            40  $^{\circ}$ C

Detector                UV (205, 235 nm)

Injection Volume      20  $\mu$ L