


# HPLC DATA SHEET

Avocado extract (아보카도 추출물, Assay)

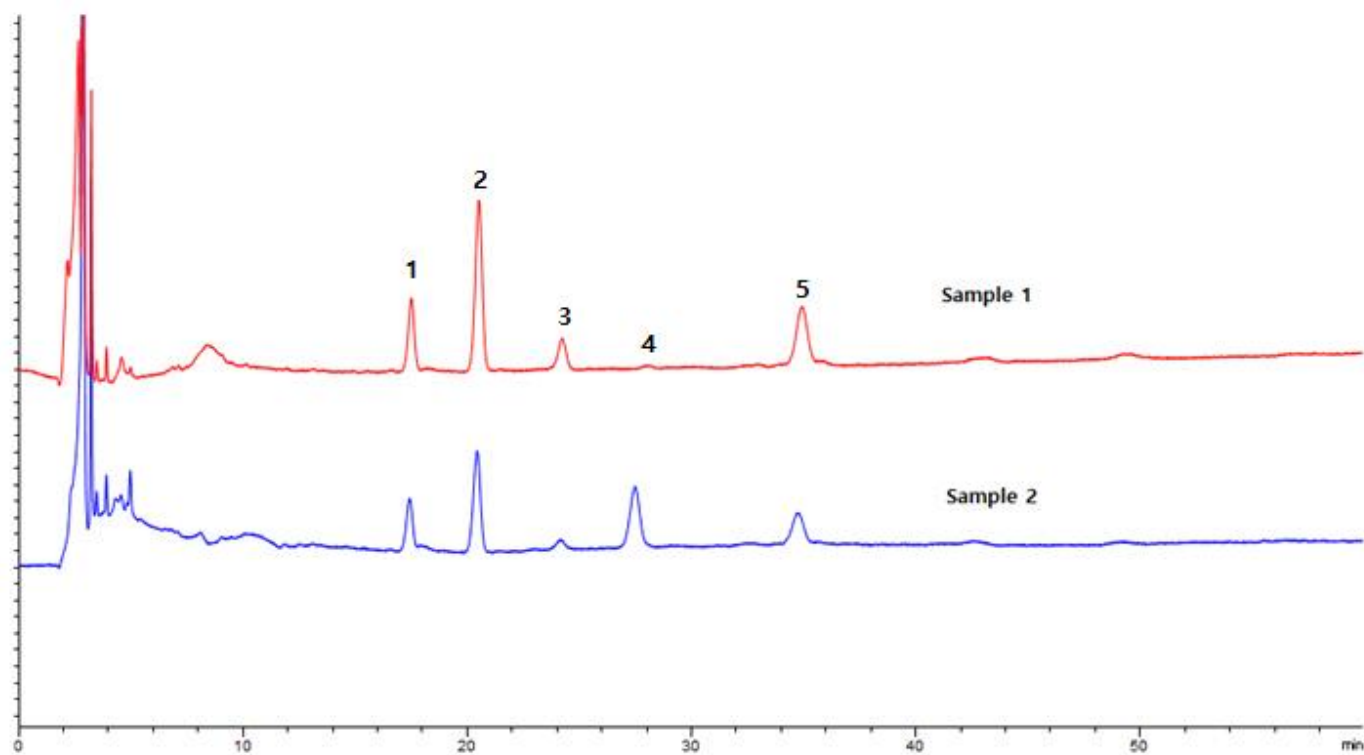

1. Ergosterol    2. Cholesterol    3. Campesterol    4. Stigmasterol    5. β-Sitosterol

## HPLC condition

Column            INNO Column C18-HC 5um, 4.6mm I.d. × 250mm (P/N 05NHC04625)  
Mobile Phase      100 % Acetonitrile  
Flow rate          1.0 mL/min  
Temperature       40 °C  
Detector           UV (278 nm)  
Injection Volume   20 uL