


# HPLC DATA SHEET

## Cholecalciferol CaCO<sub>3</sub> Tablets (Dissolution)

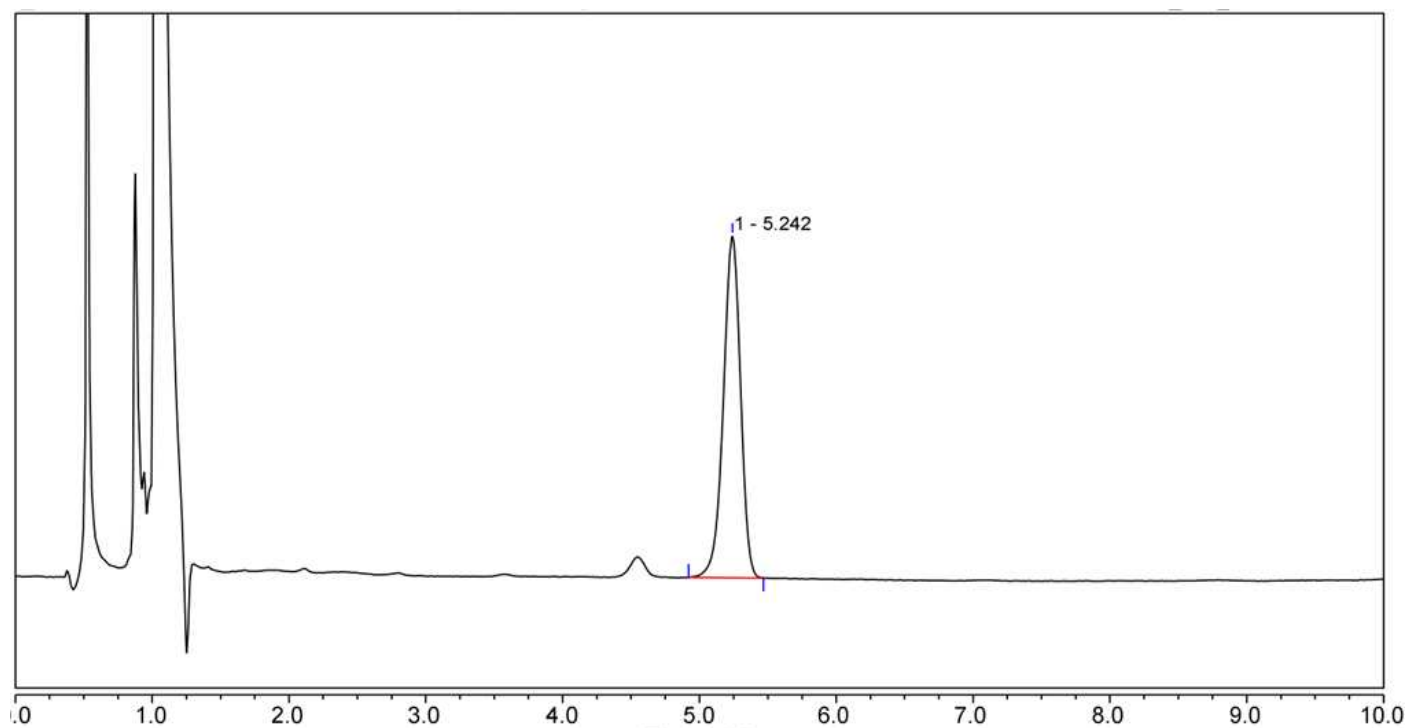

### HPLC condition

Column	INNO Column C18-NE 5 $\mu$ m, 4.6mm I.d. $\times$ 150mm (P/N 05NE04615)
Mobile Phase	Methanol : Acetonitrile (35 : 65)
Flow rate	2.0 mL/min
Temperature	30 $^{\circ}$ C
Detector	UV (254 nm)
Injection Volume	20 $\mu$ L