


HPLC DATA SHEET

Ibuprofen Capsules (KP, Related substances)

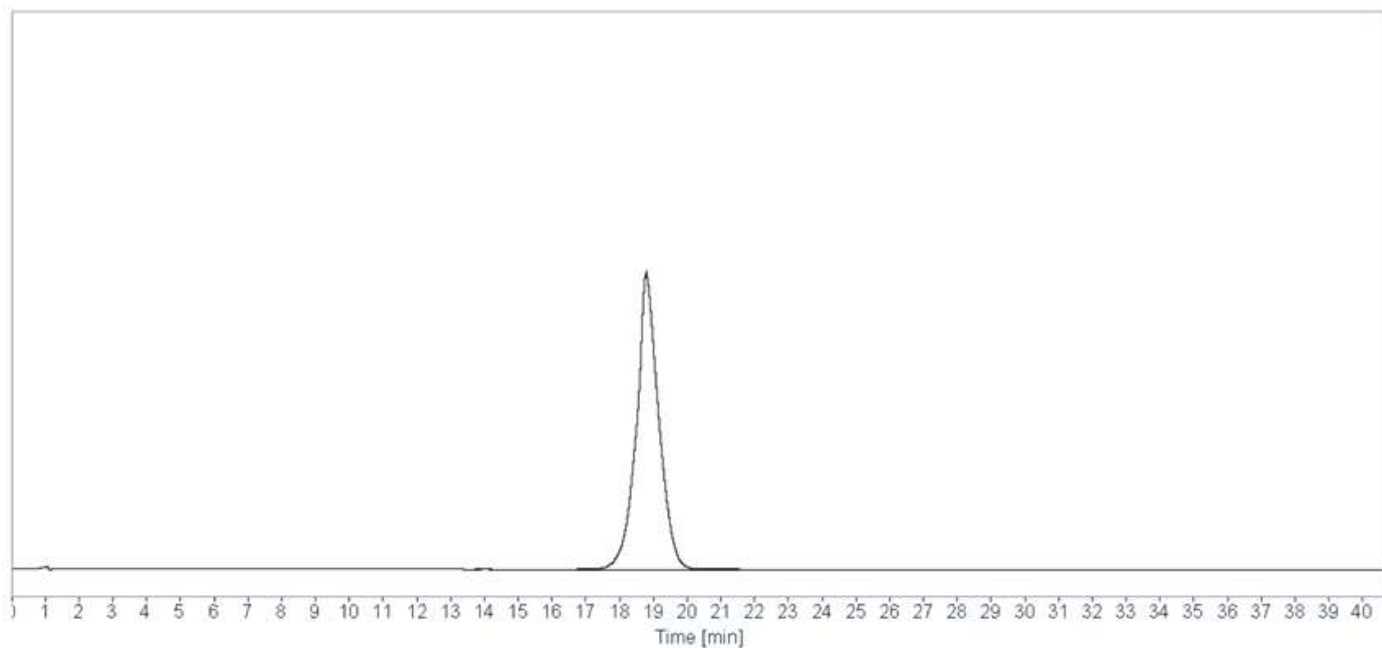

HPLC condition

Column	Aegispak C18-FA 5 μ m, 4.6mm I.d. \times 150mm (P/N 05AFA04615)
Mobile Phase	Add 0.5mL of phosphoric acid, 340mL of acetonitrile, and 600mL of water, and add water to make 1000mL
Flow rate	2.0 mL/min
Temperature	25 $^{\circ}$ C
Detector	UV (214 nm)
Injection Volume	20 μ L