


HPLC DATA SHEET

MP PP (KP, Preservative test) (1)

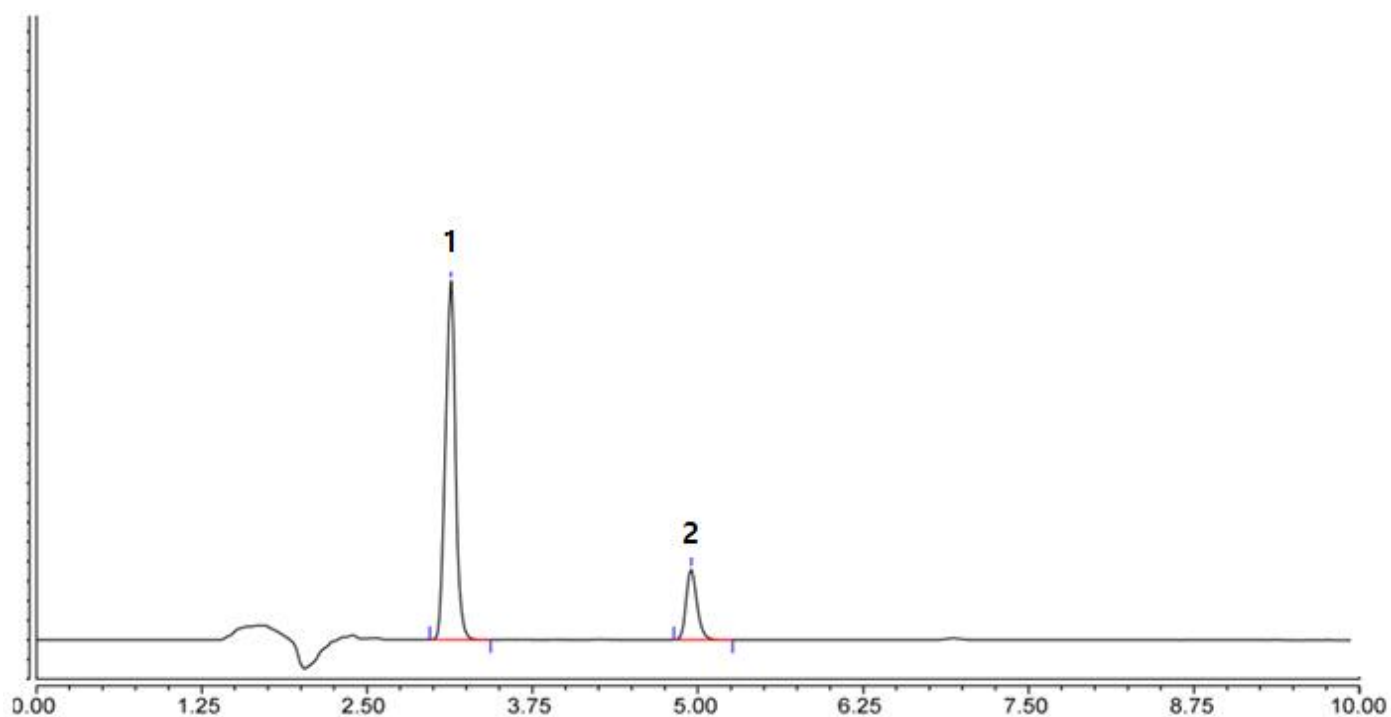

1. Methyl paraben 2. Propyl paraben

HPLC condition

Column INNO Column C18-SB 5 μ m, 4.6mm I.d. \times 250mm (P/N 05NSB04625)

Mobile Phase Acetonitrile : Water : Acetic acid anhydrous (55 : 44 : 1)

Flow rate 1.3 mL/min

Temperature 25 $^{\circ}$ C

Detector UV (254 nm)

Injection Volume 10 μ L