


HPLC DATA SHEET

Bisacodyl, Docusate Sodium Tablets (KP, Assay) (2)

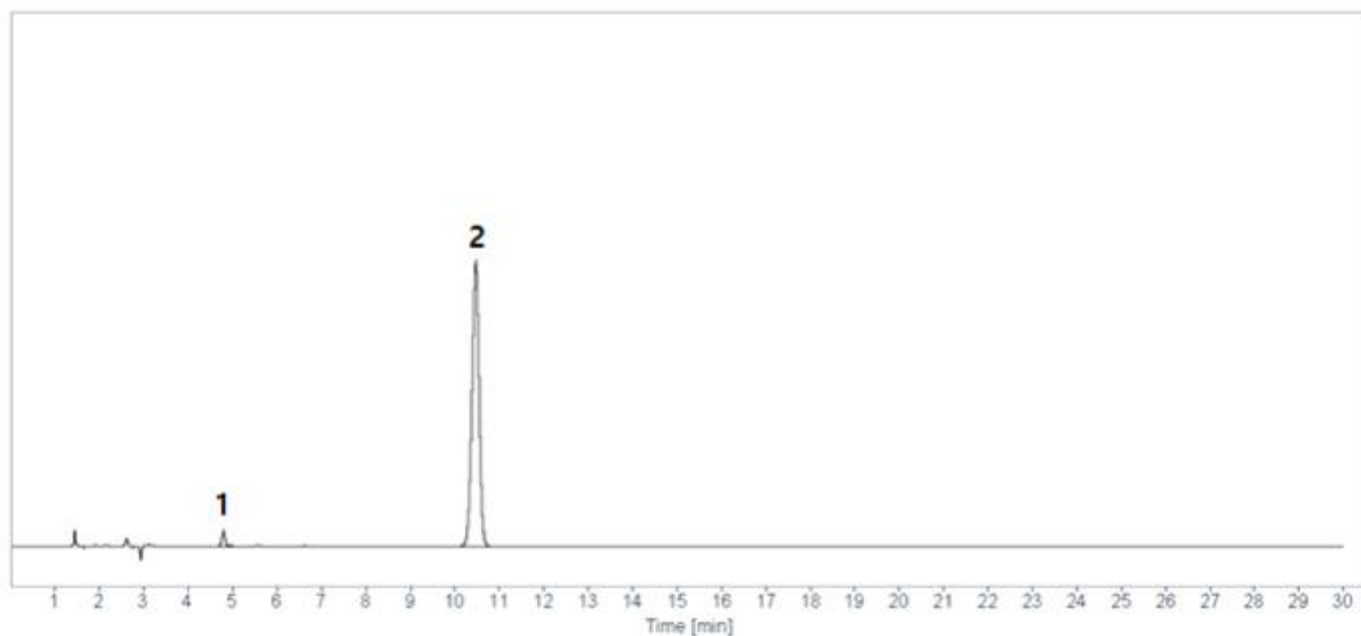

1. Docusate Sodium 2. Bisacodyl

HPLC condition

Column INNO Column C18-NE 5 μ m, 4.6mm I.d. \times 250mm (P/N 05NE04625)
0.005 mol/L Tetrabutyl ammonium phosphate : Acetonitrile : Methanol = 26 : 23 : 1

Mobile Phase - 0.005 mol/L Tetrabutylammonium phosphate : Dissolve 1.7 g of
tetrabutylammonium phosphate in 1 L of water.

Flow rate 1.0 mL/min

Temperature 25 $^{\circ}$ C

Detector UV (214 nm)

Injection Volume 10 μ L